


SolaCare - Symmetry Plus XT



MEDTRICA

Multi Zoned & Layered Foam Support

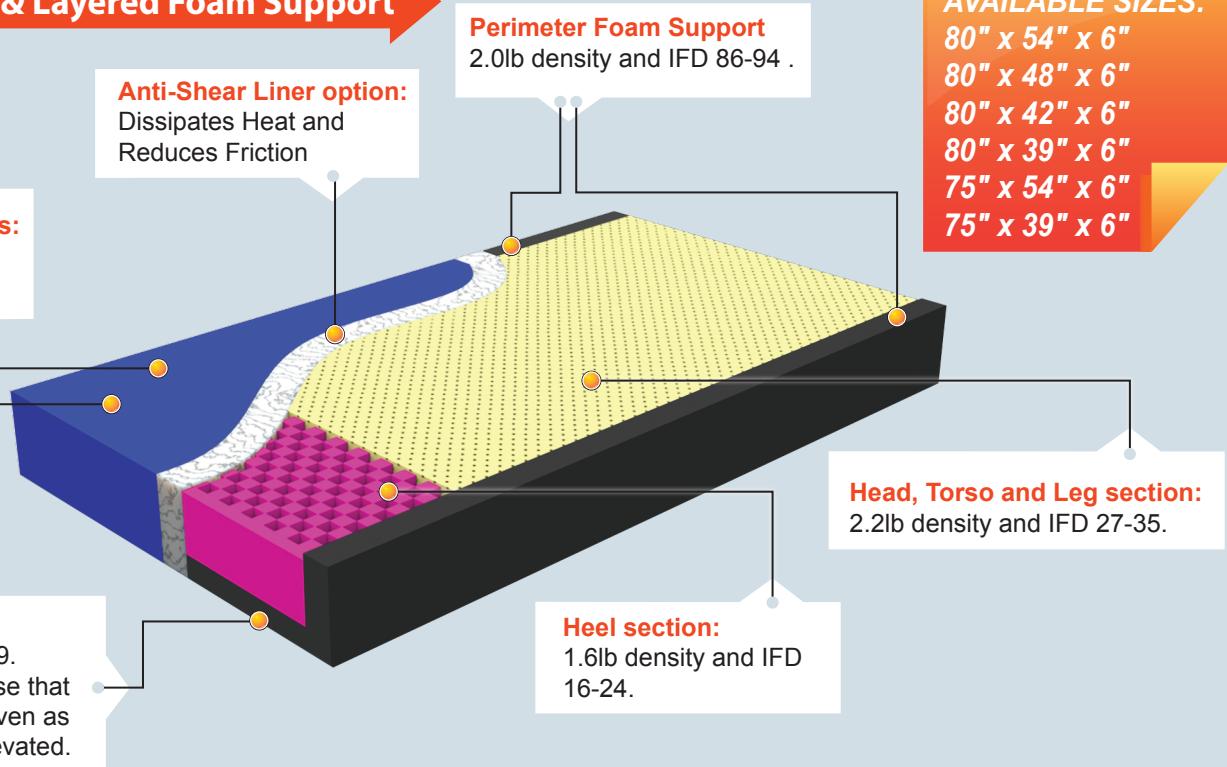

Fluid Proof Cover Options:

- Dartex™ 4 way stretch
- 70D Nylon

Anti-Shear Liner option:

Dissipates Heat and Reduces Friction

Perimeter Foam Support

2.0lb density and IFD 86-94 .

AVAILABLE SIZES:

- 80" x 54" x 6"
- 80" x 48" x 6"
- 80" x 42" x 6"
- 80" x 39" x 6"
- 75" x 54" x 6"
- 75" x 39" x 6"

Optional Sloped Heel Section

Base Foam Support

2.2lb density and IFD 51-59. Provides a firm support base that prevents "bottoming out" even as the head section is fully elevated.

Heel section:

1.6lb density and IFD 16-24.

Head, Torso and Leg section:

2.2lb density and IFD 27-35.

NON-ALLERGENIC

The entire foam core is constructed of 100% medical grade Green Friendly non-allergenic foam.

ANTI-MICROBIAL

All components are anti-microbial and assist in prohibiting growth of bacteria and fungi.

FLAMMABILITY

Meets Federal Standard 16 CFR 1632, Cal. TB 117. Optional barrier to meet Cal. TB 129 and Boston Fire IX-11.

CLAM-SHELL DESIGN COVER

Clam-shell design cover w/Double Folded Seams and 4 Reinforced Handles.

STANDARD MATTRESS COVER

Durable, snug fitting low denier 70D Nylon is soft, moisture vapor permeable, and resistant to staining with low slip vinyl bottom cover. Will not allow fluids to pass through and can be easily cleaned with standard hospital cleaning products. 100% 70D Nylon, no butyl backing.

OPTIONAL MATTRESS COVER

Dartex™ 4 way Stretchable Pressure Reducing Fabrics are the world leader in Healthcare fabrics. Dartex™ fabrics are engineered to provide the best in MVT, pressure management, abrasion and tear resistance while maintaining optimal performance throughout the entire lifespan of the product.

GUARANTEE

A full non prorated guarantee on all SolaCare Series mattresses instead of the industry standard limited warranty.

- 5 year cover guarantee
- 5 year foam core guarantee

SAFETY, COMFORT & DURABILITY

The SolaCare Symmetry Plus XT Bariatric Mattress is specifically engineered for lasting optimum performance by using better pressure reducing materials. Comfort, safety and durability is enhanced by the:
2 Layer, 2 Zone design, Firm internal foam side rails, Anti-Shear Liner and Fluid-proof Dartex™ or 70D Nylon fabric cover. Angled Heel section; No Handles, Zipper or Internal foam rail models available. The Symmetry Plus XT Bariatric accommodates patients up to 950 lbs and is Latex Free.

info@medtrica.com
www.medtrica.com



A better quality mattress, more durable, safe and comfortable - for a better price.